
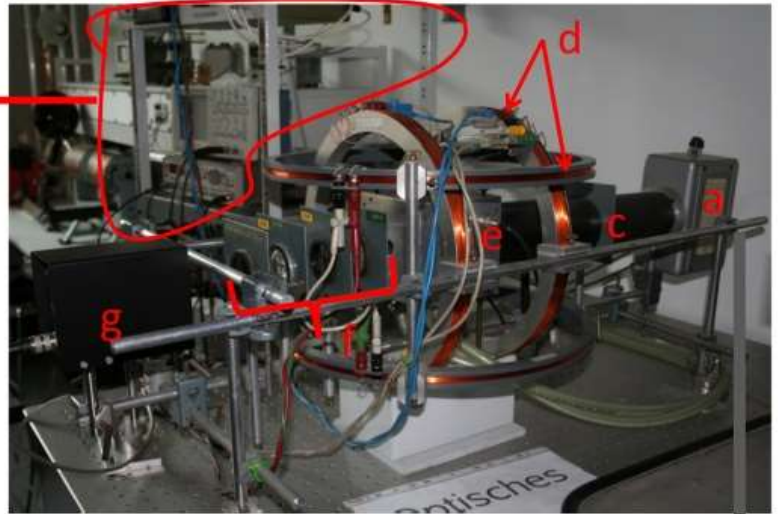


Optical-pumping Setup

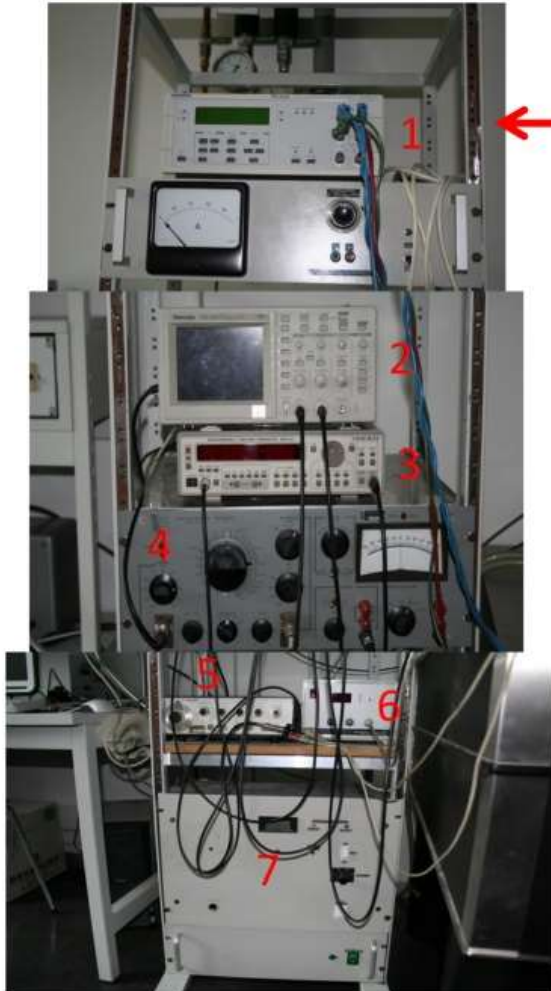

1. Current Power supply
2. Oscilloscope
3. Programmable function generator
4. Lock-in-amplifier
5. RF generator
6. Power supply of Rb lamp
7. Heating system for Rb absorption cell

- a. Detector
- b. Voltage gain for PM
- c. Focusing lens
- d. 2 pairs of Helmholtz coil
- e. Rb absorption cell
- f. Optics(Interference filter+lens+linear polarizer+ $\lambda/4$ plate)
- g. Rb lamp

Tasks:

- 1 . Investigate optical pumping of Zeeman states, measure the relaxation time and optical pumping time.

Hint: compensate earth magnetic field of the two directions ($B_v=44\mu\text{T}$, $B_h=20\mu\text{T}$).

$$B = \left(\frac{4}{5}\right)^{\frac{3}{2}} \frac{\mu_0 n I}{R}$$

n =number of windings in Helmholtz coil.

R = is radius of Helmholtz coil system.

$$\mu_0 = 4\pi \times 10^{-7} \text{ N/A}^2$$

2. Measure the resonance frequencies of Zeeman energy states, calculate Lande factor and compare with the theoretical value.

$$E = h\nu$$

$$E = g_F \left(\frac{e\hbar}{2m} \right) B m_F$$

3. Observe quadratic Zeeman Effect. (use of lock-in-amplifier)

Note (just for introduction)

Prepare the following questions:

1. what is Zeeman effect?
2. draw the S, P energy states diagram of Rb. (including Zeeman states)
3. what is the angular momentum of photon?
4. what are the electron transition selection rules between S and P states?
5. what is stimulated emission?